


Rapidpoint™ COAG




1. Swipe test card through reader

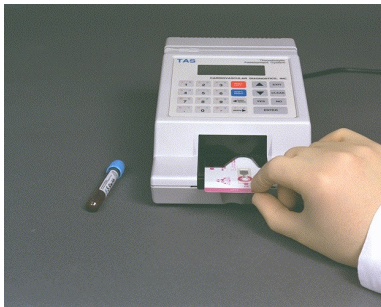

2. Insert test card in slot on analyzer

3. Select sample type , press Enter

4. Enter patient ID number, press Enter

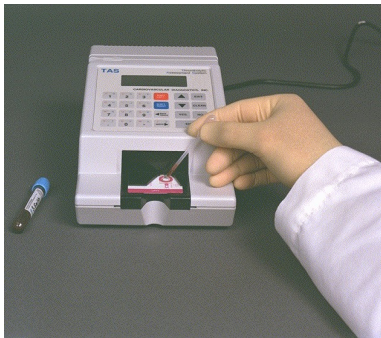

5. Add sample drop (35 uL)

6. When analyzer prompts you, remove test card

7. Press Enter at the next screen

8. Up to 1000 results are stored automatically

Daily Electronic Quality Control (EQC)

- Swipe EQC card through reader
- Place EQC card in slot on analyzer
- Analyzer will cue you to hit “Enter” after each level
- Results are stored automatically

Weekly Liquid Quality Control

- Swipe HMT card through reader
- Place HMT card in slot on analyzer
- Select “control plasma,” press Enter
- Select “level 1” or “level 2” on next screen, press Enter
- When screen says “add sample,” break QC ampoule and mix vigorously 20 times
- Discard 5 drops of liquid QC and apply 1 large drop to test card
- When analyzer prompts you, remove test card
- Results are stored on-board