


# Rapidpoint™ COAG



1. Pass card thru reader



3. Insert Card

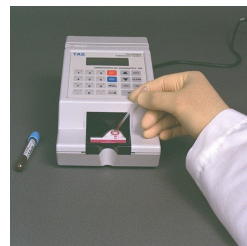

5. Enter Patient ID (ENTER)

2. Enter Your ID (ENTER)

4. Verify/Select

Sample Type (ENTER)

6. Add Sample

Bayer 

## LIQUID QUALITY CONTROL (LQC)

- EACH TEST TYPE HAS 2 LEVELS
- BREAK AMPULE AND MIX VIGOROUSLY 30 TIMES
- PREPARE AND INSERT TEST CARD FOR QC
- DISCARD 3 DROPS INTO FOIL POUCH
- HOLD VERTICALLY OVER TESTING AREA OF CARD
- APPLY ONE LARGE DROP TO SAMPLE WELL
- LOG QC RESULTS ON APPROPRIATE LOG SHEETS
- QC RESULTS STORED BY TEST TYPE

## ELECTRONIC QUALITY CONTROL (EQC)

- 3 LEVEL TEST
- EACH LEVEL HAS TIME AND SIGNAL STRENGTH

LEV	TIME	RANGE	SIGNAL
1	30	29.1- 30.9	see package
2	90	87.3 - 92.7	insert for
3	270	261.9- 278.1	range

- TO PERFORM SIMPLY SWIPE CARD THRU READER
- AND FOLLOW ANALYZER PROMPTS
- NO OPERATOR INPUT FOR EQC RANGES
- RESULTS STORED AS 3 RECORDS

[illegible]

<b>Contact Information For Additional Assistance</b>							
<b>Technical Assistance</b>					<b>1-877-BAYER11</b>		
<b>Customer Service</b>					<b>1-800-255-3232</b>		

# Bayer Diagnostics & You;

# Partners In Excellence